


Light Switch Extension Lever

USER QUICK GUIDE

Introduction

For some users with arthritis, limited grip strength or finger dexterity, it can be difficult to use a light switch. The Light Switch Extension Lever allows the user to flip the switch without needing any grip strength or nearly as much precision.

Features

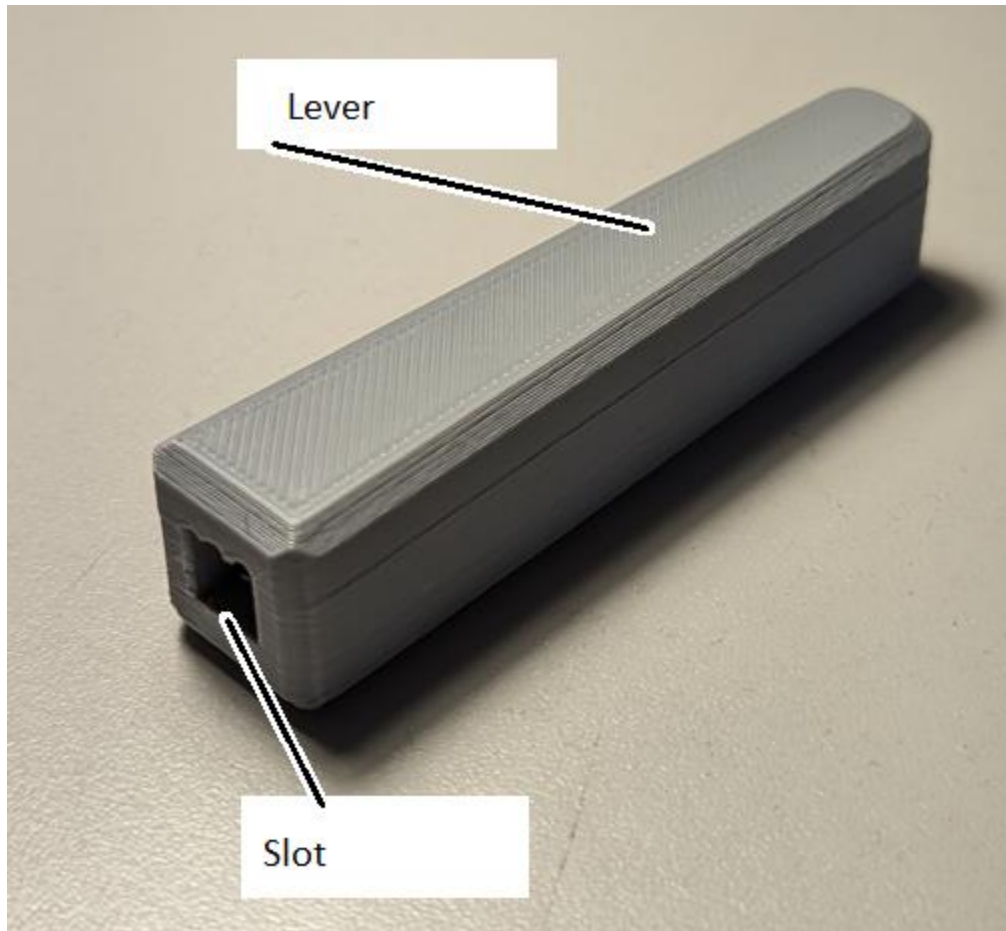

Light Switch Extension Lever

USER QUICK GUIDE

Usage

To use, place the slot over the existing toggle switch on the light switch. Push it in to the point where it wedges in place on the toggle switch. The lever can then be toggled the same way as the existing switch, but with a much larger arm needing less force and less dexterity.



Compatibility

This switch extension is compatible with any toggle style light switch and has adjustable CAD parameters to accommodate any thickness of switch toggle.

Specifications

Length: 78 mm

Thickness: 18 mm

Weight: 10 g

Cleaning

The device can be wiped clean with a damp cloth.